



*Dependable Products From People You Trust*



# **HTHP Filter Press with Threaded Cells** for Drilling Fluids Testing

**#170-181: 115 Volt**  
**#170-181-1: 230 Volt**

## **Instruction Manual**

Updated 4/14/2026  
Ver. 9

**OFI Testing Equipment, Inc.**  
11302 Steeplecrest Dr. · Houston, Texas · 77065 · U.S.A.  
Tele: 832.320.7300 · Fax: 713.880.9886 · [www.ofite.com](http://www.ofite.com)

©Copyright OFITE 2014

# ***Table of Contents***

<b>Intro.....</b>	<b>2</b>
<b>Specifications .....</b>	<b>3</b>
<b>Components .....</b>	<b>4</b>
<b>Safety.....</b>	<b>8</b>
<b>Quick Start .....</b>	<b>15</b>
<b>Operation.....</b>	<b>18</b>
<b>Data.....</b>	<b>25</b>
<b>Cell Cap Assembly .....</b>	<b>26</b>
<b>Maintenance.....</b>	<b>27</b>
<b>Diagram .....</b>	<b>30</b>
<i>Test Cell .....</i>	<i>30</i>
<b>Warranty and Return Policy .....</b>	<b>31</b>

## ***Intro***

The OFI Testing Equipment (OFITE) High Temperature High Pressure (HTHP) Filter Press is designed to evaluate the filtration characteristics of drilling fluids, cement slurries, fracturing fluids, and completion fluids under elevated temperatures and pressures. Evaluating fluids under HTHP conditions similar to the downhole environment is of paramount importance. Fluid properties must be monitored while under high temperatures and pressures as filtration behavior and wall cake building characteristics of permeable formations change with changing environments. These characteristics are affected by the shape, type, and quantities of solids present in the fluid and their physical and electro-chemical interactions, all of which are affected by changing temperatures and pressures.

OFITE manufactures and provides HTHP filtration units in two basic sizes, 175 mL and 500 mL capacities. Both are used extensively throughout the world and in all environments, but in general the 175 mL units are designed for field portability, while the larger 500 mL units are designed for laboratory usage at higher temperatures and pressures. All OFITE Filtration devices fully conform to American Petroleum Institute (API) specifications.

A complete HTHP Filter Press consists of a controlled pressure source, usually Nitrogen pressurization or Carbon Dioxide bulbs for the 175 mL units. Top and bottom pressure manifolds are provided to simulate the differential pressures found in a down-hole environment, and to prevent evaporation of the base fluid if exceeding the boiling point of that fluid. The test cells are provided in a variety of assemblies, depending upon the type of fluid tested, the filter media, and the temperatures and pressures desired. The test cells are encased inside a heating jacket, which is adjustable.

A variety of filter media are available, the most common being standard API filter paper, cement screens, and ceramic filters. The ceramic filters may be obtained to match the pore throat or permeability of the formation. Natural formation filters or cores may also be used of differing pore throat / permeability sizes. Slotted disks of varying sizes are frequently used for lost circulation materials studies.

The 175 mL has a maximum temperature of 350°F (177°C) and the 500 mL heating jackets are capable of reaching 400°F (204°C). However, lower fluid volumes due to fluid expansion at higher temperatures, limit the 175 mL units to a useful working temperature of 300°F (149°C). Anyone running tests above 350°F (177°C) must substitute a complete set of o-rings after each and every test.

## Specifications

Size:	7.5" × 11" × 23.5" (19.1 × 27.9 × 59.7 cm)
Weight:	27 lb. (12.3 kg)
Shipping Size:	20" × 13" × 13" (51 × 33 × 33 cm)
Shipping Weight:	33 lb (15 kg)
Maximum Temperature (Heating Jacket):	350°F (177°C)
Maximum Temperature (Cell):	500°F (260°C)
Maximum Pressure (Cell):	5,000 psi
Maximum Pressure (Receiver):	750 psi (5.1 MPa)
Pressure Source:	CO <sub>2</sub> Bulbs
Test Cell Capacity:	175 mL
Receiver Volume:	15 mL
Heater:	400 Watt
Power Requirement:	115 VAC, 5 Amps, 50/60 Hz 230 VAC, 3 Amps, 50/60 Hz

# Components

## Supplies:

- #153-14 Graduated Cylinder, 50 mL × 1 mL
- #154-10 Dual-Scale Thermometer with Dial, 5" Stem, 50° - 500°F (0° - 250°C)
- #165-44-2 Anti Seize Compound, Silver, 7g Pouch
- #170-19 Filter Paper, 2½" (6.35 cm), Specially Hardened for Filter Presses, 100/Box
- #170-35 Adjustable Wrench, 6"

## Assemblies:

### #170-04 CO<sub>2</sub> Pressurize Unit:

- #143-02-10 CO<sub>2</sub> Puncture Head Assembly
  - #143-02-12 Puncture Pin
  - #143-02-13 O-ring
  - #143-02-14 O-ring
- #143-03 Barrel for CO<sub>2</sub> Cartridge
- #170-08-004 Regulator
- #170-20 Manifold Block
- #170-32 Needle Valve, Male, ⅛" × ⅛"
- #171-23-1 Safety Pin with Lanyard
- #171-34 Gauge, 1,500 psi, 2" Face, ¼" NPT Bottom

### #170-06 Back Pressure Receiver, 15-mL Stainless Steel Tube for CO<sub>2</sub>

- #143-00-001 Regulator
- #143-01 Gauge, 200 psi, ⅛" Bottom Connection
- #143-02-10 CO<sub>2</sub> Puncture Head Assembly
  - #143-02-12 Puncture Pin
  - #143-02-13 O-ring
  - #143-02-14 O-ring
- #143-03 Barrel for CO<sub>2</sub> Cartridge
- #143-06 Safety Bleeder Valve
- #144-11 Elbow
- #170-07 O-ring
- #170-28 Receiver Body
- #170-32 Needle Valve, Male, ⅛" × ⅛"
- #171-23-1 Safety Pin with Lanyard

**#170-181-S Test Cell, Stainless Steel (For Mud Testing)**

#120-910-028	O-ring for Rupture Disk, Viton 75D, Qty: 1
#130-81-040	Retaining Ring, Qty: 2
#165-44-2	Anti Seize Compound, Silver, 7g Pouch, Qty: 4
#170-13-3	O-ring for Cell, Viton 75D, Qty: 4
#170-16	Valve Stem, Qty: 2
#170-17	O-ring for Valve Stem, Viton 75D, Qty: 4
#170-18	Cement Screen, Qty: 2
#170-180-020-S	Cell Body, 175 mL, Qty: 1
#171-190-023	Locking Ring, Qty: 2
#171-190-027	Rupture Disk, Qty: 1
#171-190-029	Cap Wrench, Qty: 1
#171-190-032-S	Cell Cap, Outlet, Cement, Qty: 1
#171-190-033-S	Cell Cap, Inlet, Cement, Qty: 1
#171-190-057	O-ring for Valve Stem, Viton 90D, Qty: 4
#171-190-058	O-ring for Rupture Disk, Viton 90D, Qty: 1
#171-190-060	O-ring for Cell, Viton 90D, Qty: 4

**#170-00-1 Heating Jacket, 115 Volt:**

**#170-01-1 Heating Jacket, 230 Volt:**

#164-32	Male Connector for Power Cable (230 Volt)
#170-05-001	Thermostat
#170-10	Thermostat Pilot Light
#170-11	Heating Element, 200W, Qty: 2
#170-15	Base
#170-21	Stand
#170-30	Stainless Steel Thermostat Cover
#170-44	Rubber Foot ½", Qty: 4
#171-32	Midget Knob
#171-82	Power Cord with Male Plug, 8' (115 Volt)

**Optional:**

#143-05	CO <sub>2</sub> Bulbs, 8-Gram, Package of 10, UN 1013
#152-00	Hamilton Beach Mixer, With Container
#155-20	Timer, 60 Min. Interval
#170-03	Carrying Case, Stainless Steel
#170-13	O-ring for Test Cell, NBR/Nitrile (Buna N), For temperatures up to 250°F (121°C)
#170-13-4	O-ring for Test Cell, Perfluorocarbon (FFKM), For temperatures up to 500°F (260°C)
#170-13-5	O-ring for Test Cell, Ethylene propylene (EPM/EPDM), For temperatures up to 400°F (204.4°C), Water-based fluids only
#170-33	HTHP Cell Cap Puller
#170-40	Test Cell Removal and Carrying Tool
#170-91	HTHP Pressure Relief Tool
#170-92	Safety Clamp for HTHP Fluid Loss Cells

**#170-181-SP Spare Parts Kit:**

Part Number	Description	Qty.
#140-60-01	O-ring for Bleeder Valve	2
#143-01	Gauge, 200 psi, 1/8" Bottom Connection	1
#143-02-13	O-ring for Puncture Pin Holder	2
#143-02-14	O-ring for Puncture Pin Holder, Rear	2
#143-05	CO <sub>2</sub> Bulbs, 8 Gram, Package of 10, UN 1013	60
#153-14	Graduated Cylinder, 50 mL × 1 mL	2
#154-10	Dial Thermometer, 5" Stem, 50° - 500°F and 0 - 250°C	1
#170-13-3	O-ring for Cell, Viton 75D	50
#170-16	Valve Stem	4
#170-17	Valve Stem O-ring	100
#170-19	Filter Paper, 2½" (6.35 cm), Hardened for Filter Press	5
#171-190-057	O-ring for Valve Stem, Viton 90D	100
#171-190-060	O-ring for Test Cell, Viton®	50
#171-23-1	Safety Pin with Lanyard	1

**Note**

Spare parts listings are intended to be used as a reference for future purchases. Everyone's consumable requirements will be different, and replacement quantities needed will depend upon the number of test performed on a daily and/or weekly basis.

**Optional Items for HTHP Filtration Testing:**

The items listed below are not included in the HTHP Filter Press, but they are items that will enable the technician to perform a more uniform and reproducible test while maintaining a high degree of safety. As optional items, the usage is not compulsory, but consideration should be given to these items when running tests at elevated temperatures and pressures.



Interval Timer, 60 minute  
(#155-20)



Cell Cap Removal Tool  
(#170-33)  
(Set Screw Cell Assemblies Only)



Cell Carrying Tool (#170-40)



HTHP Pressure Relief Tool  
(#170-91)  
(To release trapped pressure)



Stand for HTHP Cell Assembly  
(#171-190-028)



Safety Shield (#171-06)

# Safety

## Pressure and Temperature Considerations

Do not use or reconfigure this equipment in a manner not specified in this manual.

### Pressure

There are two reasons for operating at elevated pressures when performing a filtration analysis.

1. To test fluids at temperatures above the boiling point, the vessel must be pressurized, which in turn elevates the vapor pressure (boiling point) so that the fluid remains liquid and does not turn to steam.
2. If testing for drilling applications, pressurization will more approximate down-hole conditions, enabling the technician to match both bottom hole temperatures and pressures.

HTHP Filter Presses are pressurized either with Carbon Dioxide or with Nitrogen gas. The OFITE 175 mL HTHP Filter Press, 170-00 Series, is small enough to be portable, so it is usually pressurized with small (65 mm length) Carbon Dioxide (CO<sub>2</sub>) bulbs, which contain 10 cm<sup>3</sup> of CO<sub>2</sub> gas and weigh 8 grams. These bulbs are pressurized to approximately 1,000 psi and they contain plenty of carrier gas to run a complete 30 minute filtration test, if running the standard API 500 psi differential test, which is usually 600 psi on the top manifold and 100 psi on the bottom or back pressure.

For temperatures less than 200°F (93°C) a Back Pressure Receiver is not required as the filtrate will not reach the boiling point. However when operating above the boiling point of water, a suitable Back Pressure Receiver is required, otherwise the test fluid will turn to steam and the test is invalidated. The standard receiver tube supplied with the 175 mL HTHP Filter Press holds 15 mL of filtrate, so it is very important that the receiver outlet ball valve is opened after approximately ten seconds of filtrate collection, as a fluid with poor filtration qualities can easily fill the 15 mL receiver tube almost instantaneously. If this happens the filtrate hopefully will all be ejected from the safety bleeder valve, but if there is too much fluid volume, the liquid may end up inside the regulator rendering it useless which will require servicing by a knowledgeable technician.

### Temperature

Normally when one experiences a loss of pressure it is not due to a failure of the metal alloy in the cell, but rather is a failure of the o-ring or elastomer, which provides the seal. These o-ring may deform or melt under elevated temperatures causing a pressurization failure of the cell, which is often sudden and catastrophic. For example, if the valve stem o-ring suddenly fails, then steam at 400°F and under great pressure may shoot outward horizontally in one or several directions. A safety shield should always be used when operating any HTHP Filter Press and especially when one is going to extreme temperatures and pressures. Above 350°F (176°C) All o-rings must be replaced after each and every test.

## Heating Jacket

Turn the thermostat knob on the front of the metal plate to change the temperature of the heating jacket. The thermostat cover is marked from 1 to 10 and each whole number represents a separate temperature. Each filter press will heat up a little differently so it is a good idea to mark with a pen directly on the face-plate of the thermostat cover the temperature which corresponds with each number. The actual temperature of the heating jacket is measured with a stem thermometer and then the thermometer is placed into the cell body.

The heating jacket will get very hot during the test so be careful not to touch the outside of the jacket at any time during the heat-up period or the test. It is especially easy to get burned when installing or removing the Back Pressure Receiver as the technician is working in a very confined space. The power cord is supplied for either 115 Volts or 230 Volts and due to the various types of plug connectors around the world it may be necessary to change the plug on the end of the power cord. The power cable is 8 feet in length and the heating jacket should be placed no farther than 8 feet from the appropriate electrical outlet. The heat up time will vary from one instrument to the next and the pilot light will turn on when the set temperature has been reached. The temperature of the fluid inside the cell however will not be at the set temperature, so always allow one hour of heating time for the fluid, after the cell has been fully inserted into the heating jacket.



The HTHP Filter Press is a very safe instrument to operate as long as the instructions are strictly followed. It is imperative that the technician thoroughly understands the assembly and disassembly procedures and also why the order of steps is in place. Care should be taken to never use defective parts and that OFITE's temperature, pressure, and fluid volume limits are strictly adhered to at all times.

The heating jacket and cell will get very hot during the normal duration of the test so the operator should take care to avoid burns. The operator should be careful when attaching and removing pressure manifolds from the heating jacket, especially the Back Pressure Receiver assembly, which involves operating in a small confined space.

### Safety Retainer Pin

Always use a Safety Retainer Pin with an attached lanyard, and get in the habit of securing the cotter pin to the pin. Improper placement of retainer pins in the past have resulted in the pressure manifold separating from the filter press which can cause damage to the equipment.



Retainer Pin with Lanyard



171-23-1 Safety Retainer Pin with Lanyard

Normally the 175 mL HTHP Filter Press uses CO<sub>2</sub> bulbs, for both portability and pressurization. Even though they are small and appear innocuous, they can contain pressures up to 1,000 psi (6,900 kPa). CO<sub>2</sub> bulbs should never be allowed to overheat. Exposure to fire or carelessly leaving a box of bulbs on an automobile dashboard in the sun may result in an explosion. Nitrous Oxide (NO<sub>2</sub>) which is available in the same sized bulbs as the CO<sub>2</sub>, should never be used as a pressure source for HTHP filtration. Under high temperatures and pressures, Nitrous Oxide can detonate in the presence of grease, oil or carbonaceous material.



CO<sub>2</sub> Bulbs, 10/Box, (#143-05)

### Working Pressure

OFITE clearly stamps the working pressure of the cell assembly along with the assembly serial number on each cell body and cell cap. These pressure limits should never be exceeded on the pressure manifolds and fluid volumes inside the cell should strictly comply with instructions (see page 19). However, in the past HTHP cells and caps from various manufacturers were often not stamped, or in some cases, very old cell assemblies were stamped "2500". Cell assemblies with set screw secured caps should never be taken to 2,500 psi, under any temperature or condition.

### O-ring Seals

Two sets of o-rings are included with the Filter Press. The first set is made from Viton 75D. These o-rings are black and should be used for tests up to 400°F only. The second set is made from Viton 90D. OFITE suggests that all o-rings (cell caps, cell bodies, and valve stem) are changed each time after running a test of 350°F (176° C) or above.

### Fluid Expansion Space

Cells should never be filled completely with fluid. Always leave some room for the fluid to expand under temperature. This keeps the cell caps from becoming overstressed and possibly leading to a cell pressure failure or even a loss of the cell cap completely.

Periodically check all gauges to ensure the proper pressures are in the cell, as the pressure often increases or decreases depending upon the test procedure and conditions.

API Recommended Void Space		
Fluid / Temperature	Void Space	Fluid Volume
Water-based drilling fluid. < 300°F	0.6" (1.5 cm)	138 mL
Water-based drilling fluid. > 300°F	1.5" (4.0 cm)	86 mL
Oil-based drilling fluid. < 350°F (176°C)	1" (2.5 cm)	115 mL
Oil-based drilling fluid. > 350°F (176°C)	Not Recommended- Use 500 mL Cell	



Void Space in a 1,500 psi rated cell body

### **Cell Corrosion**

Test fluids under high temperature and pressure can corrode the cell body and caps. Carefully inspect the cell body and caps for corrosion before and after each test.


Some materials are more susceptible to corrosion than others. Also, some fluids and additives are more corrosive than others. OFITE offers a variety of cell materials for different levels of corrosion resistance and cost.

### **Regulator Safety**

The OFITE 175 mL HTHP Filter Press uses two regulators as standard equipment. Both may be operated with Carbon Dioxide (CO<sub>2</sub>) or Nitrogen (N<sub>2</sub>) gas.

1. Top or Main Pressure:
  - a. Maximum Inlet Pressure - 3,000 psi
  - b. ¼" NPT Connections
  - c. Single stage regulator
2. Bottom or Back Pressure
  - a. Maximum inlet pressure - 2,500 psi
  - b. ¼" NPT Connections
  - c. Single stage regulator
3. Never subject the regulator to inlet pressure greater than its rated inlet pressure as shown on the regulator body.
4. Do not allow oil, oil bearing materials, grease, or other combustibles to contaminate the inside or outside of the regulator, especially where Oxygen is in use.
5. To minimize heat effects, avoid the use of piping or tubing between the pressurized cylinder and the regulator. If unavoidable it should be as short as practical and equipped with a shutoff valve located just before the regulator.
6. It is not recommended to pressurize HTHP Filtration tests with Oxygen due to its potentially explosive nature.
7. Never pressurize a regulator that has loose or damaged parts or is in questionable condition. Never loosen a connection or attempt to remove a part until gas pressure has been relieved.
8. Before transporting pressurized cylinders that are not secured on a cart designed for such transport, remove all regulators and recap the cylinders.



9. Check the regulator and all connections for leaks after installation and periodically thereafter. Also check for leaks after any service in which parts or connections were loosened. Brush with an approved leak detection solution and bubbles should indicate a leakage.
10. An appropriately sized pressure relief device downstream of the regulator should be installed in your system to prevent damage to equipment and/or injury to personnel should an internal failure of the regulator occur.
11. Do not purge oxidizing or flammable gases in the presence of flame, lit cigarettes, or other sources of ignition or towards people.

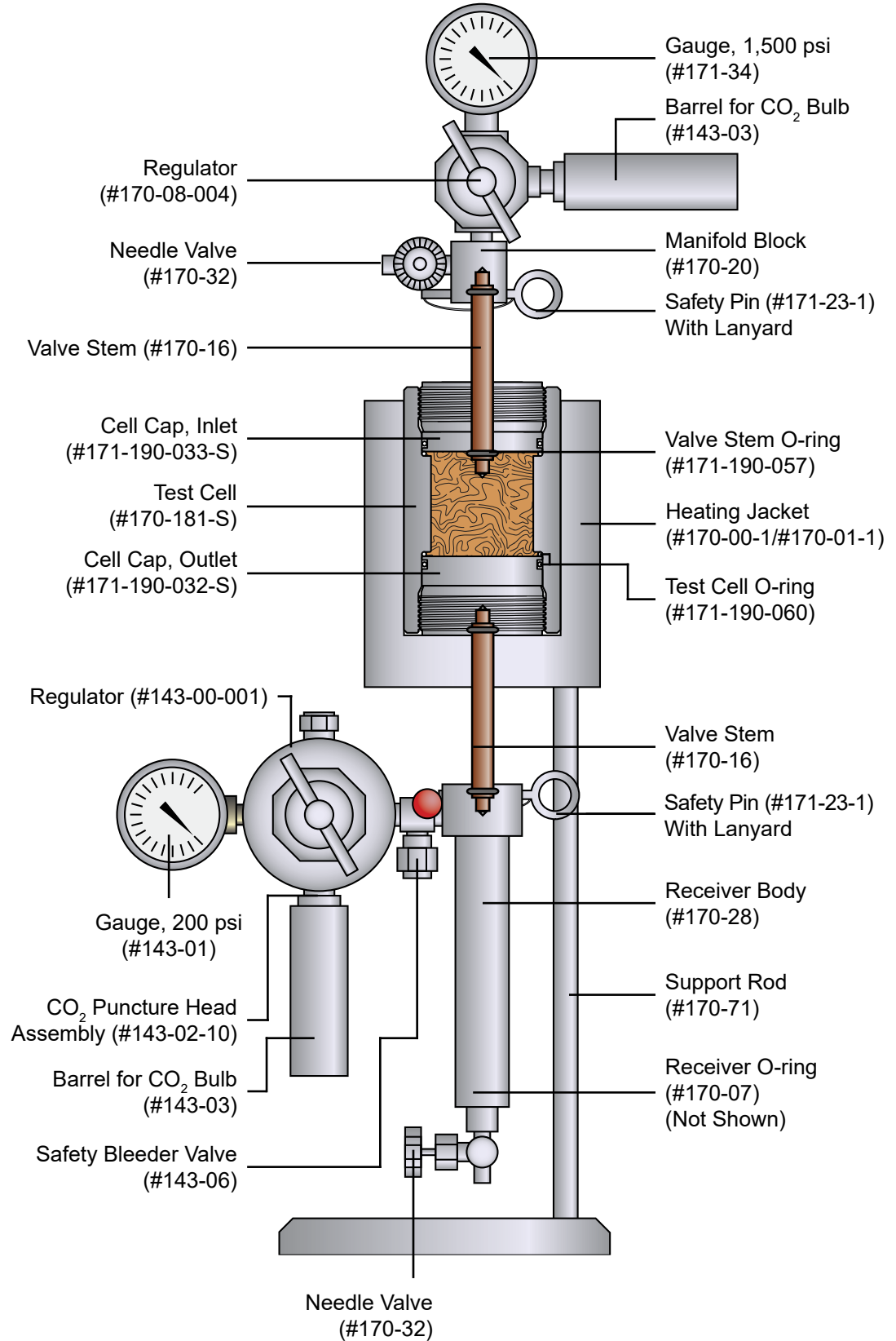

## Quick Start

1. Place the thermometer in the heating jacket.
2. Preheat the heating jacket to 10°F (6°C) above the test temperature. The pilot light will turn on at the target temperature.
3. Inspect all o-rings. Replace any that show signs of damage or wear.
4. Apply grease to all o-rings.
5. Assemble both cell caps. Place an o-ring in the groove around the cap.
6. Make sure all threads are clean and free of debris.
7. Apply anti seize compound (#165-44-2) to the threads of the cell locking rings and valve stems.
8. Screw the inlet cell cap into the cell and tighten it completely.
9. Screw a valve stem into the inlet cell cap and tighten it completely.
10. Pour the test fluid into the cell. Leave adequate space for fluid expansion.
11. Place an o-ring in the groove inside the cell.
12. Place a circle of filter paper on top of the o-ring.
13. Place an o-ring in the groove around the cell cap.
14. Screw the outlet cell cap into the cell. Tighten it completely.
15. Screw a valve stem into the outlet cell cap and tighten it completely.
16. Invert and place the cell inside the heating jacket and rotate the cell until it is seated. Start the hour heat-up time.
17. Transfer the thermometer to the cell body.
18. Connect the Inlet (top) pressure assembly and lock it in place with the retaining pin.
19. Connect the back pressure assembly and lock it in place with the retaining pin.
20. Puncture CO<sub>2</sub> bulbs and adjust the inlet/back pressure regulators to a desired back pressure.
21. Open (loosen) the inlet (top) valve stem ½ turn to pressurize the cell.
22. Heat for one hour. The thermometer should register a target temperature. The total heat up not to exceed 1 hour.

23. Increase pressure on top pressure unit to 500 psi over the back pressure.
24. Open bottom valve stem to initiate filtration.
25. Collect and note filtrate amount at 10 seconds, 1.0, 7.5, and 30.0 minutes.
  - a. 10 second collection is cautionary as a fluid with poor filtration control can quickly fill up the receiver.
  - b. If back pressure raises above 100 psi reduce the pressure by opening the ball valve.
26. Collect filtrate for a total of 30 minutes, maintaining desired inlet and back pressures.
27. After 30 minutes close (tighten) inlet and outlet valve stems to seal off the cell.
28. Turn top and bottom regulator T-screws counterclockwise to release pressure.
29. Open the outlet valve on the back pressure receiver to collect all the filtrate in the graduate cylinder.
30. Open needle or bleeder valves on both the top and back pressure units.
31. Remove the top and bottom retainer pins and remove pressuring units from cell.
32. Drain any residual filtrate collected from the top of the receiver into the graduate cylinder.
33. Allow the cell to cool in the heating jacket.
34. Once the cell has cooled to room temperature, remove it from the heating jacket.
35. Open (Loosen) the inlet valve stem to relieve pressure on the cell.
36. Unscrew the outlet cell cap.
37. Wash the filter cake and record the thickness to the nearest  $\frac{1}{32}$ ".
38. Clean and dry the filter press jacket, cell, and pressure assemblies thoroughly.
39. Inspect and replace any or all o-rings (All if test run > 300°F (149°C)).
40. Report: Double the filtrate volume - Correct to the API standard of 7.1 in<sup>2</sup>. Spurt Loss volume is Optional.

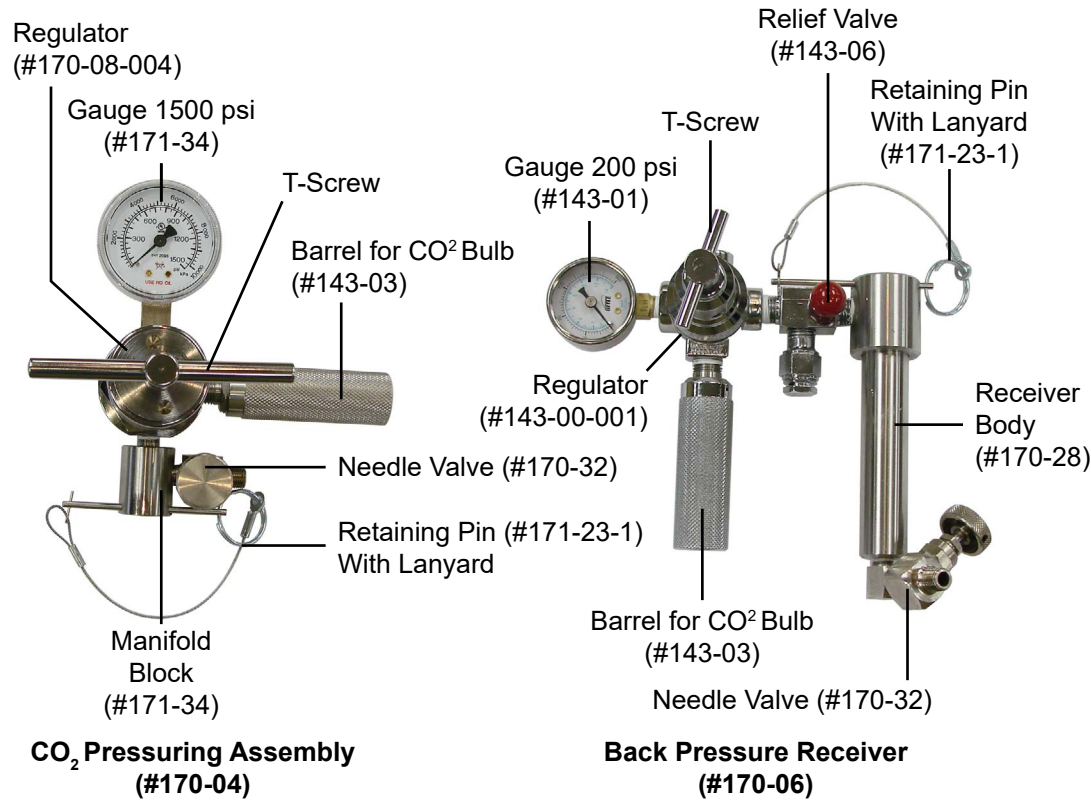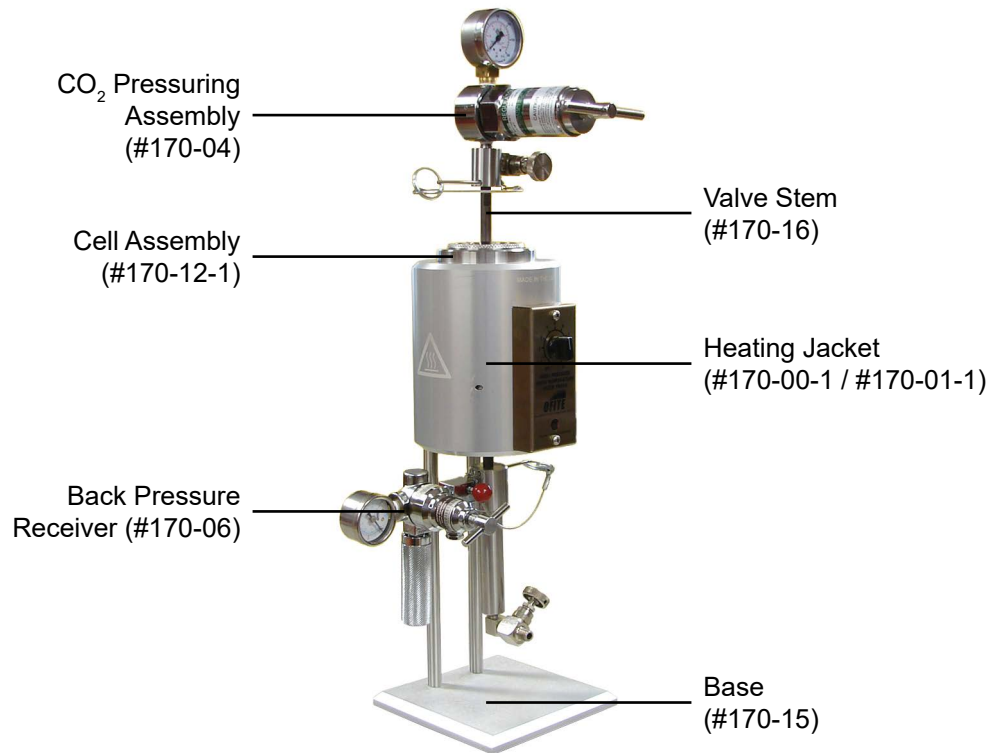





- Screw the inlet cell cap into the cell body.

The cell cap should turn smoothly in the threads and not require any tools for complete tightening. If you encounter resistance, carefully unscrew the cap and start over. Make sure the threads are properly seated before tightening completely.

- Place an o-ring in each of the two o-ring grooves on both valve stems.
- Screw one valve stem into the inlet cell cap. Tighten the valve stem completely.



To increase the life of the valve stem and cap, apply a thin layer of high-temperature thread lubricant (#165-44-1) to the tip of the valve stem.

- Prepare the test fluid.
- Invert the cell body and carefully pour the sample into the cell. Leave enough void space for thermal expansion (see chart below).

API Recommended Void Space		
Fluid / Temperature	Void Space	Fluid Volume
Water-based drilling fluid. < 300°F	0.6" (1.5 cm)	405 mL
Water-based drilling fluid. > 300°F	1.5" (4.0 cm)	353 mL
Oil-Based Drilling Fluid. Any temperature.	1" (2.5 cm)	382 mL



Be careful not to spill fluid on the o-ring shoulder inside the cell.

- Place an o-ring in the cell and another on the cell cap.
- Place a circle of filter paper or ceramic disk on top of the cell o-ring. Gently push the paper or disk downward so it contacts the o-ring without the paper binding or pinching.
- Screw the outlet cap into the cell body.
- Screw the other valve stem into the outlet cell cap and tighten it completely.
- Invert the cell and place it in the heating jacket with the outlet (filter) side pointed down. Rotate the cell in the heating jacket so that the pin in the bottom of the heating well seats into the hole in the bottom of the test cell. This will anchor the cell in the well and prevent it from rotating as the valve stems are opened and closed.




**Note**

18. Move the thermometer from the heating jacket to the thermocouple hole on top of the test cell. The heating time of the sample should be one hour (no more and no less). Set a timer to 60 minutes. Begin the heating at this time.

Do not initiate filtration when the thermometer indicates the desired temperature has been reached. The fluid inside the cell may not have yet reached the proper temperature. An API funded project found that for most fluids, it usually takes upwards of an hour for the fluid inside the cell to reach the target temperature. If accurate temperature measurement is required during the heating portion of the test, a Thermocouple Assembly (#171-45-1) is available which will directly and accurately measure the fluid temperature inside the cell.

19. While the cell is heating, connect the CO<sub>2</sub> pressuring assembly to the inlet valve stem and lock it in place with the retaining pin and lanyard. Check that the safety bleeder valve is closed (knob tightened clockwise).
20. Connect the back pressure receiver to the bottom valve stem and lock it in place with the retaining pin and lanyard. Ensure that the safety bleeder valve is closed (knob pushed all the way inward) and the receiver outlet needle valve is closed (knob tightened clockwise). Be careful. It is easy to get burned in the confined space below the heating jacket.

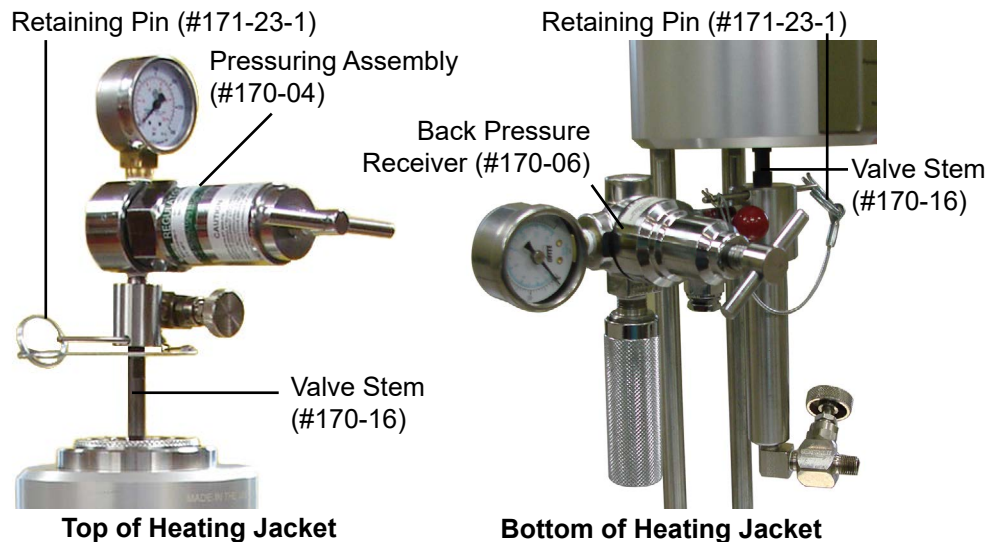

21. Keeping the valve stems closed, adjust the top and bottom regulators to the recommended back pressure for your test.



**Note**

The upper and lower limits of the test pressure differential are determined by the test temperature. As this temperature exceeds 212°F (100°C), the back pressure must be increased in order to prevent vaporization of the filtrate. If running a standard API test, the 500 psi differential pressure must be maintained, so the top pressure will have to be increased accordingly. The table below shows the pressures recommended for various test temperatures.

<b>Recommended Minimum Back Pressure</b>					
Test Temperature		Vapor Pressure		Minimum Back Pressure	
°F	°C	psi	kPa	psi	kPa
200-299	95-149	12.1-67	84-462	100	700
300-374	150-189	67-184	462-1269	200	1400
375-399	190-199	184-247	1269-1704	275	1900
400-424	200-219	247-326	1704-2245	350	2500
425-450	220-230	326-422	2245-2912	450	3100

22. Open (loosen) the Inlet valve stem ½ turn to pressurize the sample. Maintain this pressure on the fluid until the temperature has stabilized. The heating time of the sample should be one hour. This will prevent vaporization of the test solution if heating above the boiling point of the liquid.
23. After the one hour heat up time, increase the pressure on the inlet pressure unit to 500 psi (3,448 kPa) more than the back pressure.
24. Open (loosen) the bottom valve stem ½ turn to initiate filtration as soon as possible.



**Tip**

Closely monitor the pressure gauges. If at any time during the test the pressure inside the cell rises above the setpoint, carefully open the needle valve on the top pressure assembly just enough to bleed off the excess pressure. Then close the valve. If the cell pressure decreases due to collection of filtrate, increase the pressure with the inlet regulator.

25. To collect filtrate, carefully open the needle valve on the bottom of the back pressure receiver while holding a graduated cylinder to the valve port. Close the valve immediately after the outlet pressure just begins to decrease and filtrate is collected.

Collect filtrate at 10 seconds, 1 minute, 7.5 minute, and 30 minute intervals. The initial 10 second collection is precautionary in nature, as a fluid with little filtration properties may fill up the 15 mL receiver tube almost immediately, potentially damaging the regulator. Do not record the 10 second volume collected, but do record the volume collected at the other time intervals. Timers pre-set to 7.5 minutes and 30 minutes are helpful.



**Note**

While collecting filtrate, maintain the test temperature within  $\pm 5^{\circ}\text{F}$  ( $\pm 3^{\circ}\text{C}$ ). If the back pressure rises during the test, cautiously reduce the pressure by opening the valve on the receiver and drawing off some of the filtrate into the graduated cylinder.

26. After 30 minutes the test is over. Immediately close (tighten) the top and bottom valve stems to seal off the cell.
27. Allow a minimum of 5 minutes for the filtrate in the receiver to cool to avoid vaporizing. Collect filtrate by opening the needle valve on the receiver and leave it open until fluid ceases to evacuate the valve.
28. Close (tighten) the needle valve on the receiver.
29. Turn the regulator T-screws on both the inlet and back pressure receiver counterclockwise until the T-screw feels "loose". This will close off the flow of pressurized gas.
30. Release the pressure from the Inlet pressuring unit to the cell by opening the needle valve (turn counterclockwise).
31. Release the pressure from the back (bottom) pressuring unit to the receiver body by pulling outward on the knob which will open the bleeder valve.
32. Remove the top and bottom retaining pins and remove the two pressure assemblies.



Take care to avoid touching the hot heating jacket, especially when removing the back pressure receiver.

33. Open the outlet needle valve and drain any residual filtrate collected in the receiver into the graduated cylinder.
34. Record the total amount of filtrate collected in the graduate cylinder by reading from the bottom of the meniscus curve. Multiply this amount in milliliters by 2.
35. Allow the cell to cool to room temperature.
36. Remove the cell from the heating jacket after once again checking that the valve stems are tightly closed. An optional Cell Carrying Tool (#170-40) makes this a simple and safe operation.



DO NOT remove the heated cell using pliers or an open ended wrench attached to the valve stem. This is very dangerous and could cause serious damage.



Important

The test cell will still be under approximately 600 psi (4.140 kPa) pressure. To avoid possible serious injury, position the cell with the outlet end down and cool it to room temperature before disassembling. The cell must be below 100°F (38°C). Cool for at least one hour at room temperature before removing the caps.

37. Place the cooled cell upright with the outlet side down. The optional stand for the HTHP cell assembly (#171-190-028) makes a great cell holder. This will help preserve the filter cake. Open (loosen) the Inlet valve stem to bleed off pressure from the cell body. There should be a noticeable release of pressure from the cell.



Stand for HTHP Cell Assembly (#171-190-028)

Hold a rag or piece of paper over the valve stem in order to catch any liquid or solids that might be ejected under pressure from the cell. This will help prevent an unsightly residue build up on the ceiling which may occur over time. It may be necessary to return the cell to the heating jacket and use the stop-pin in order to loosen the valve stem.



Note

Pressure should NOT be relieved from the cell by opening the outlet valve stem as the filter cake may seal off the cell.

Safe Release  
of Pressure  
from a Cell



A good indication that pressure is trapped inside the cell body occurs if the cell caps are unusually difficult to loosen. To relieve any trapped pressure, an OFITE Pressure Relief Tool (#170-91) may be inserted through the inlet cell cap hole which will relieve any trapped pressure inside the cell. Ensure this function is performed on the inlet side of the cell as the filter cake may seal off the outlet end of the cell.



Important

38. Unscrew and remove the cell caps.

39. There are several ways to remove the filter along with the filter cake.
  - a. Hold the cell and invert it over a sink and let the weight of the fluid inside the cell push the disk out.
  - b. Careful manipulation with a spatula will often remove a filter and filter cake intact without losing fluid.
40. Discard the fluid inside the cell unless it is required for further testing. Save the filter cake for analysis.
41. Clean and dry the apparatus thoroughly after each use.
42. Inspect all o-rings for deformities and replace as necessary. Replace all o-rings if the test was conducted above 350°F (177°C).



Tip

After each test, the HTHP filter press should be left so that a new test may be performed with no clean-up and a minimum of assembly.

## Data

### Filtrate Volume

The HTHP filter press has a filtration area of 3.55 in<sup>2</sup> (22.9 cm<sup>2</sup>). This is half the area of a standard filtration test, which is 7.1 in<sup>2</sup> (45.8 cm<sup>2</sup>). To compare the results of this test to a standard filtration test, double the total filtrate volume collected.

$$V_F = 2 (V_{30})$$

Where:

$V_F$  = Standard Filtrate Volume (mL)

$V_{7.5}$  = Filtrate volume collected after 7.5 minutes

### Spurt Loss (Optional):

Spurt Loss is the amount of filtrate collected before the filter cake has had a chance to form and is expressed in millimeters. To calculate the spurt loss, use the following equation:

$$V_1 = 2 [V_{7.5} - (V_{30} - V_{7.5})] = 2 (2V_{7.5} - V_{30}) = 4V_{7.5} - 2V_{30}$$

Where:

$V_1$  = Spurt Loss

$V_{7.5}$  = Filtrate volume collected after 7.5 minutes

$V_{30}$  = Filtrate volume collected after 30 minutes

### Filter Cake

Wash the filter cake on the paper with a gentle stream of water. Measure and report the thickness of the filter cake to the nearest 1/32 in (0.8 mm). A ruler with the "zero mark" at the very edge of the ruler is useful here. Cake descriptions may be subjective and such notations such as hard, soft, rubbery, and fine, etc. convey adequate information on cake quality.

# Cell Cap Assembly



Tip

1. Choose the appropriate cap for your test:
  - 171-190-030-S - Outlet, with 60 Mesh Screen for Filter Paper
  - 171-190-031-S - Inlet, with 60 Mesh Screen
  - 171-190-034-S - Outlet, Scribed for Ceramic Disks
2. Place the locking ring (#171-190-023) around the cap.
3. Place the retaining ring (#130-81-040) into the groove around the outside of the cap. Make sure it engages completely around the circle.
 

The cap should turn freely inside the locking ring.
4. Place an o-ring in the port in the cap. Wrap a rupture disk (#171-190-027) with nickel anti-seize tape (#171-190-040) and screw it into the port.




Locking Ring  
(#171-190-023)

Retaining Ring  
(#130-81-040)

Rupture Disk  
(#171-190-027)

Cap

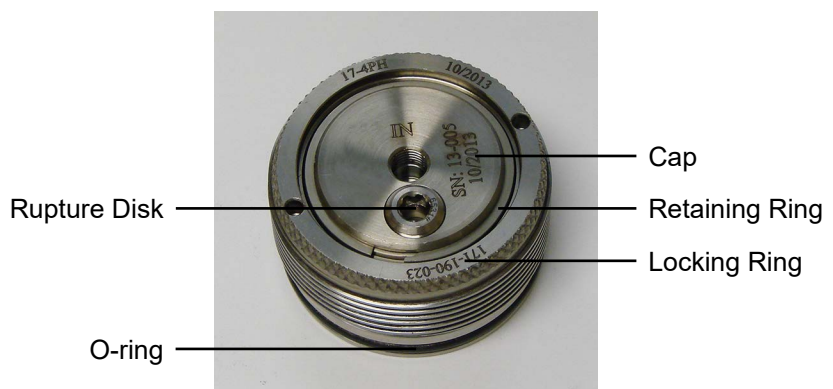


# Maintenance

1. HTHP cell assemblies should periodically be checked for:
  - a. Cell Corrosion – Stress cracking or pitting
  - b. End Cap Compression
  - c. Elastomers – Check for flat spotting and melting
  - d. Valve stems should be sharply pointed and not dull on outlet side.
2. Always thoroughly clean up after each tests.
  - a. Replace all elastomers if previous test was in excess of 300°F (149°C).
  - b. Clean out and blow dry if possible, the back pressure receiver.
  - c. Cell body must be completely cleaned with special attention to the o-ring groove.
3. Periodically check for leaks.
  - a. Cell Assembly – Immerse sealed pressurized cell assembly in a sink filled with water – look for bubbles.
  - b. Top Pressure Assembly – Immerse regulator, valves, and sealed cell assembly in a sink full of water and look for bubbles. Hold assembly by the gauge and do not let the gauge get wet.
  - c. Back Pressure Receiver – Similar to “b”. Hold assembly by the gauge and do not allow it to get wet.

4. Regulator and Manifold Maintenance:

- a. Replace both o-rings on the puncture pin (#143-02-13 and #143-02-14) at least once a year.

**#142-10 - CO<sub>2</sub> Pressuring Assembly**

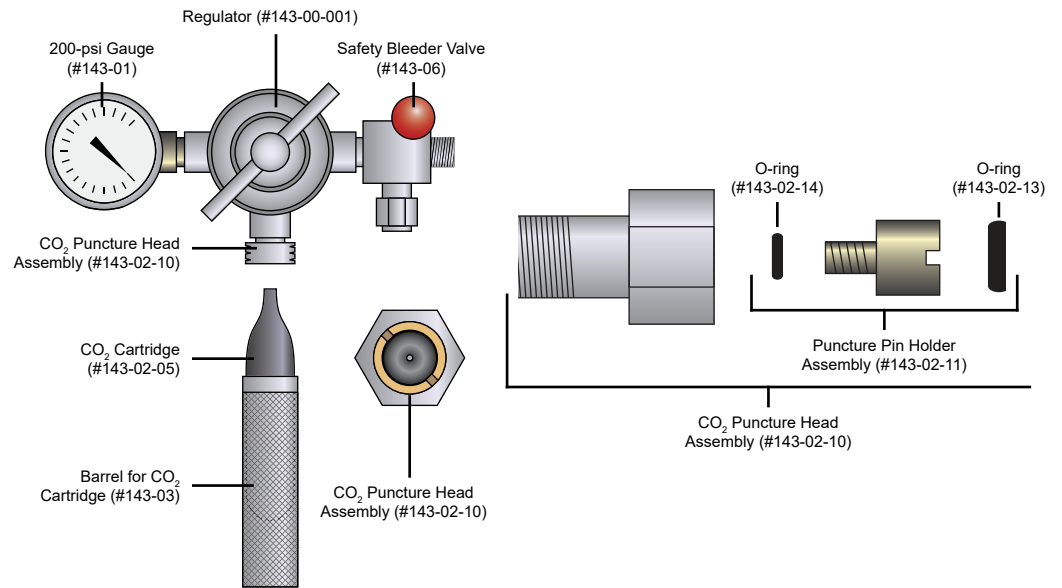

- b. If the regulator loses pressure or steadily increases in pressure, repair the regulator.

5. Power Cord

- a. Check power cord for insulation wear and loose connections near the heating jacket and plug.
- b. The heating jacket should be placed no farther than the distance of the power cord from an electrical outlet.
- c. The power cord should be kept away from the hot surface of the heating jacket, while in use.




Any customer installed wiring, power cords or electrical connectors will void all warranties.

## Regulator Maintenance and Trouble Shooting

Symptoms	Cause: Resolution
Gas leak at the regulator outlet when the adjusting screw is loosened fully counterclockwise	Seat leak or *creep: Repair the regulator
Outlet pressure increases while downstream valves are closed	Seat leak or *creep: Repair the regulator
Gas leak from the spring housing case	Diaphragm failure: Repair the regulator
Excess drop in outlet pressure with the regulator flow open	Blocked seat assembly or inlet filter: Repair the regulator
Gas leak from any pipe thread joint	Loose fitting: Remove the connection. Clean the affected surfaces. Reapply Teflon tape and tighten.
Gas leak from relief valve	Faulty relief valve: replace the valve. Seat leak or *creep: Repair the regulator
Inconsistent repeat readings	Seat sticking: Repair the regulator. Bad pressure gauge: Replace the gauge.
Gauge does not return to zero with no pressure applied to the regulator.	Gauge has physical damage: Replace the gauge.
<p>*Creep is an increase in outlet pressure that occurs when pressure escapes even when the valve is closed. Regulator seats can be compromised by particles in the process stream which can cause minor imperfections in the sealing surface. The high flow and small orifice created during pressure regulation combine to turn a very small particle into a fast projectile. This projectile can nick the sealing surface of the seat and cause leaks. Filtering particulates from the process stream should be a high priority, and a small filter can reduce the potential for creep and increase the life expectancy and accuracy of the regulator.</p>	

# Diagram

Test Cell

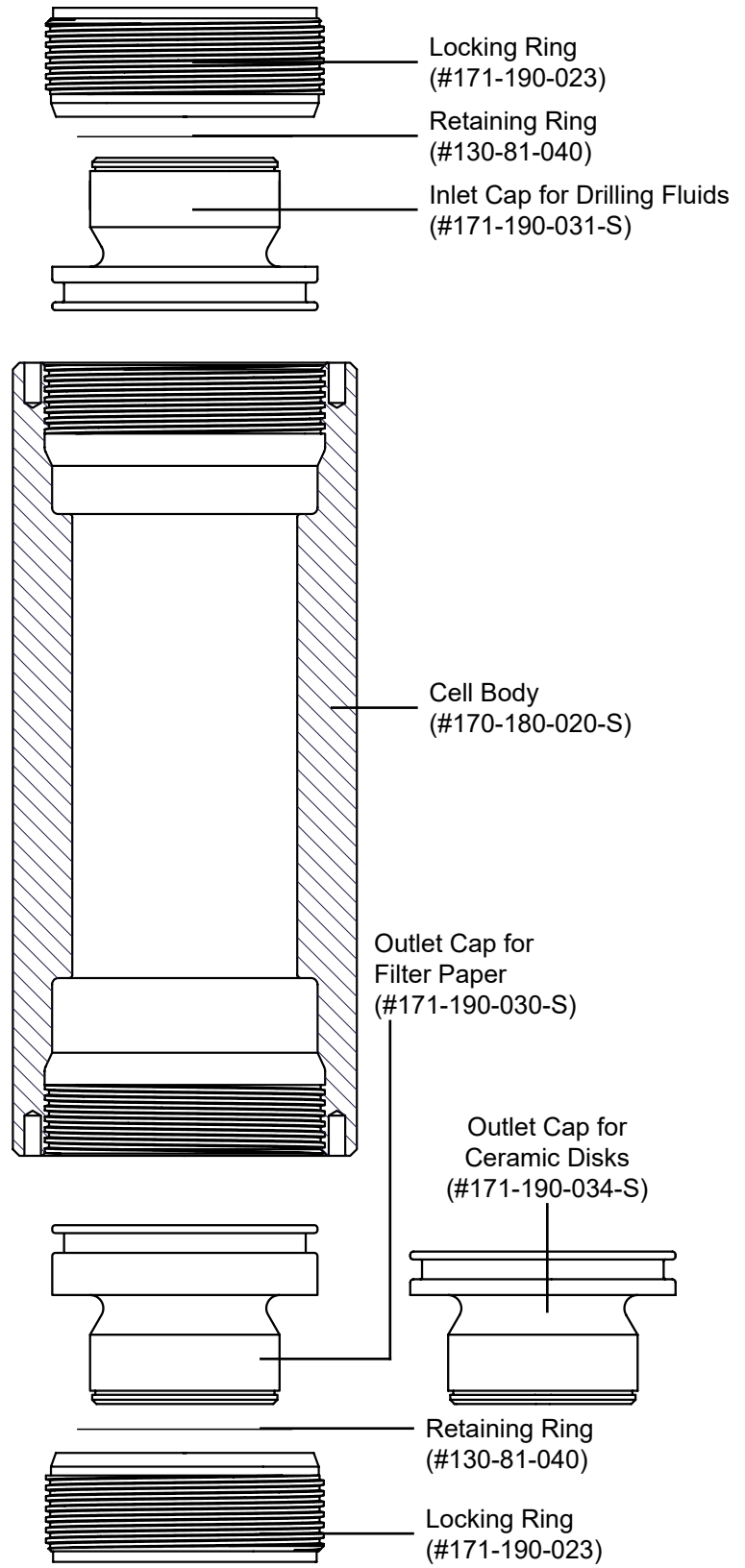

# Warranty and Return Policy

## Warranty:

OFI Testing Equipment, Inc. (OFITE) warrants that the products shall be free from liens and defects in title, and shall conform in all respects to the terms of the sales order and the specifications applicable to the products. All products shall be furnished subject to OFITE's standard manufacturing variations and practices. Unless the warranty period is otherwise extended in writing, the following warranty shall apply: if, at any time prior to twelve (12) months from the date of invoice, the products, or any part thereof, do not conform to these warranties or to the specifications applicable thereto, and OFITE is so notified in writing upon discovery, OFITE shall promptly repair or replace the defective products. Notwithstanding the foregoing, OFITE's warranty obligations shall not extend to any use by the buyer of the products in conditions more severe than OFITE's recommendations, nor to any defects which were visually observable by the buyer but which are not promptly brought to OFITE's attention.

In the event that the buyer has purchased installation and commissioning services on applicable products, the above warranty shall extend for an additional period of twelve (12) months from the date of the original warranty expiration for such products.

In the event that OFITE is requested to provide customized research and development for the buyer, OFITE shall use its best efforts but makes no guarantees to the buyer that any products will be provided.

OFITE makes no other warranties or guarantees to the buyer, either express or implied, and the warranties provided in this clause shall be exclusive of any other warranties including ANY IMPLIED OR STATUTORY WARRANTIES OF FITNESS FOR PURPOSE, MERCHANTABILITY, AND OTHER STATUTORY REMEDIES WHICH ARE WAIVED.

This limited warranty does not cover any losses or damages that occur as a result of:

- Improper installation or maintenance of the products
- Misuse
- Neglect
- Adjustment by non-authorized sources
- Improper environment
- Excessive or inadequate heating or air conditioning or electrical power failures, surges, or other irregularities
- Equipment, products, or material not manufactured by OFITE
- Firmware or hardware that have been modified or altered by a third party
- Consumable parts (bearings, accessories, etc.)

## Returns and Repairs:

Items being returned must be carefully packaged to prevent damage in shipment and insured against possible damage or loss. OFITE will not be responsible for equipment damaged due to insufficient packaging.

Any non-defective items returned to OFITE within ninety (90) days of invoice are subject to a 15% restocking fee. Items returned must be received by OFITE in original condition for it to be accepted. Reagents and special order items will not be accepted for return or refund.

OFITE employs experienced personnel to service and repair equipment manufactured by us, as well as other companies. To help expedite the repair process, please include a repair form with all equipment sent to OFITE for repair. Be sure to include your name, company name, phone number, email address, detailed description of work to be done, purchase order number, and a shipping address for returning the equipment. All repairs performed as "repair as needed" are subject to the ninety (90) day limited warranty. All "Certified Repairs" are subject to the twelve (12) month limited warranty.

Returns and potential warranty repairs require a Return Material Authorization (RMA) number. An RMA form is available from your sales or service representative.

Please ship all equipment (with the RMA number for returns or warranty repairs) to the following address:

OFI Testing Equipment, Inc.  
Attn: Repair Department  
11302 Steeplecrest Dr.  
Houston, TX 77065  
USA

OFITE also offers competitive service contracts for repairing and/or maintaining your lab equipment, including equipment from other manufacturers. For more information about our technical support and repair services, please contact [techservice@ofite.com](mailto:techservice@ofite.com).